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CE

1. Introduction

We would like to thank you for choosing a MOBIAK SA product. The product user manual is an important document that provides you with information on how to use it and technical features. **Please read the user manual carefully before using the product. Learn the product and practice before any application.**

- This manual applies to our model: **0805905.**
- The annotations and illustrations in this manual might be slightly different with the real parts due to quality improvement or changing design. Please check it carefully.
- The product you have in your hands may differ from the one in the following photos. The assembly process is similar.

SYMBOLS: Please refer on ANNEX A

2. Intended Use:

The 0805905 bed is used for patients who have health problems and need to stay in bed. It is designed for people who need to be hospitalized for health or disability reasons in nursing homes, hospitals or healthcare facilities. It was also designed as a comfortable solution for home care of people who need care. Its purpose is to support this care.

3. CAUTION (!):

- In case of damage or malfunction of your product DOT NOT use the product and please contact the authorized and trained representative.
- Do not use the product for a purpose not indicated in this manual.
- MOBIAK S.A declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration to the frame of the product.
- MOBIAK S.A reserves the right to change the information contained in this document without previous notice.
- There is a ± 3%. Deviation for the Dimensions of the Products.

4. GENERAL WARNINGS:

- If you don't read this user manual, it is preferable not to use this product or another available part.
- If you don't understand the usage or the precautions, please contact with the dealer or the proper technical person before you use the product because may be caused damage.
- Keep away the product from heat sources.
- Save this document for future reference.
- Be careful when children are nearby and do not allow the children to play with the product.
- Do not exceed the maximum load capacity.
- Special care is always required where there are moving parts that can cause limb entrapment and injury.
- Do not attempt to lift the product by any removable parts.

5. Product Presentation:

- Electric Hospital Bed.
- Ideal for Intensive Care Units.
- Helps the patient benefit and have better care during recovery from cardiac, respiratory, and other procedures.
- With one button, it enables caregivers to convert the hospital bed into a chair.
- Integrated bed controller on both side rails, convenient for long-term operation for both the patient and the caregiver.
- Both electric and manual CPR, provide quick release to a flat position for emergency cardiopulmonary resuscitation.
- Electric leg section extension.
- Angle indication, both for back, Trendelenburg & Anti-Trendelenburg.
- Fifth wheel for direction, anti-vibration for better transportation.
- Central locking on both sides of the bed.
- With emergency battery backup.
- With night light, on the control unit.
- With stainless steel push handle
- With four bumpers on the four corners.
- Electric backrest lift control

- Electric leg lift control
- Electric height adjustment control
- Electric trendelenburg control
- Electric anti-Trendelenburg control
- Rails with easy height folding mechanism, for immediate access to the patient
- With 2 pcs of 1.3Ah lead-acid battery, with an autonomy of approximately 5-7 full* cycles of movement operation
- The time needed to fully charge the battery is approximately 12 hours.
- The battery should be fully recharged every 3 months, in order to maintain its capacity
- In the event of a power outage, the battery will automatically power the system. The standard activation procedure is then followed.
- Low-vibration back recline function to reduce forces exerted on the abdomen
- Auto contour function to check the correct patient position
- Removable head and footboard with bumpers to prevent impact and protect against impacts - four bumpers at the four corners.
- Castors with central brake, rotating
- With fifth wheel, anti-vibration for direction and better transport.
- Wide height adjustment range from 56cm to 98cm
- Wired controller
- Built-in Function Control on the external part of the siderails for use by healthcare personnel
- Built-in Patient Function Control (in the internal part of the siderails)-(control functions can be locked via the built-in control for use by healthcare staff)
- Manual release for CPR on both sides
- CPR button on the control panel
- Sheet-bed linen base (located at the bottom-base of the bed and can be used by nursing personnel when changing sheets)
- Oxygen cylinder base suitable for 2,3 liters cylinders
- Integrated slots for placing drainage bags, urine collectors
- Night light for better patient service during the night
- Both electric and manual CPR provide quick release to a flat position for emergency cardiopulmonary resuscitation.
- Central locking on both sides of the bed
- With 2 hole bases for placing I.V stands
- With 2 hole bases for placing lifting poles
- With an adjustable I.V stand, with 4 folding hooks
- With lifting pole

**Full operating cycle: The bed performs all of its movements, back and leg lifting, height adjustment, Trendelenburg & Reverse Trendelenburg, bed-to-chair conversion function.*

6. Mattress Selection

The mattress must be purchased separately. Use mattresses or air mattresses according to the opinion of the attending physician. The size of the mattress should not be larger than the bed frame. The mattress **0805941 HOSPITAL MATTRESS "PRO I"** with dimensions 191 x 82 x 14cm is suitable for the size of the bed.

7. Specifications

Model	0805905
Description	HOSPITAL BED
Functions	Electrical & Manual
Frame	Steel
Bed platform	Steel
Overall Length (without bed extension)	202cm
Overall Length (with bed extension)	222cm
Overall width (with siderails up)	100cm
Overall width (with siderails down)	98cm
Mattress platform	196 x 90cm
Minimum Height	56 cm
Maximum Height	98 cm
Backrest raising tilt	0-75°
Legrest raising tilt	0-45°
Trendelenburg	0-15°
Reverse Trendelenburg	0-15°
Castors	4 pcs , 15cm , with brakes
Fifth wheel	1 pc, 10cm
Max load	250 Kg
Voltage	100-240V
Frequency	50/60 Hz
Hospital grade wire and plug Water resistance degree	IPX4
Product Net weight	250kg
Frame painting	Electrostatic
Motors	5


8. Precautions For Use


- Do not exceed the maximum load capacity for any reason.
- When lowering the head or foot section, make sure that nobody inserts hands, legs or other body parts under the bed base.
- Make sure that the floor surface supporting the bed is flat**.**


9. Before Every Use:


- Check if the frame is damaged in order to guarantee a safe use of the product. (There are no cracks or fractures in the frame).
- In case of damage do not use the product and contact your dealer for further instructions.
- Check, if the product is properly assembled and all screws are secure and well screwed on.
- Always check the condition of wear for the mechanical parts to ensure the absolute safe use of the product for people and objects.
- In case of damage do not use the product and contact your dealer for further instructions.

Pictures of Correct and Incorrect Use of Product

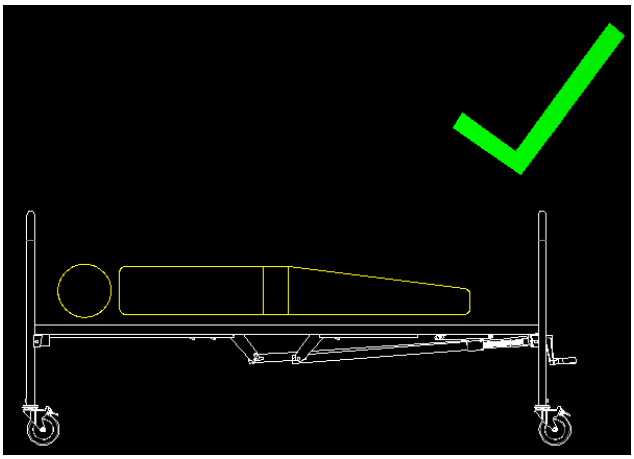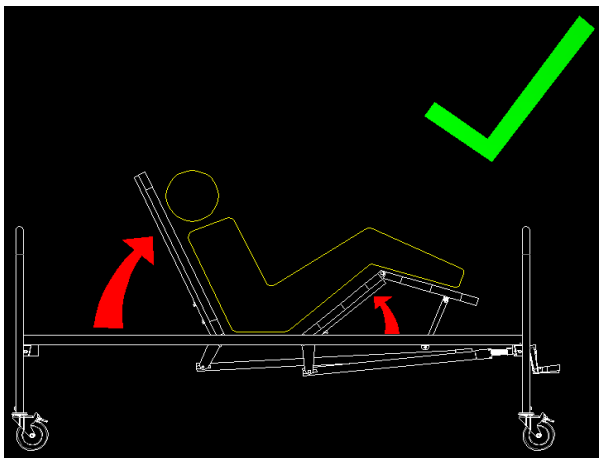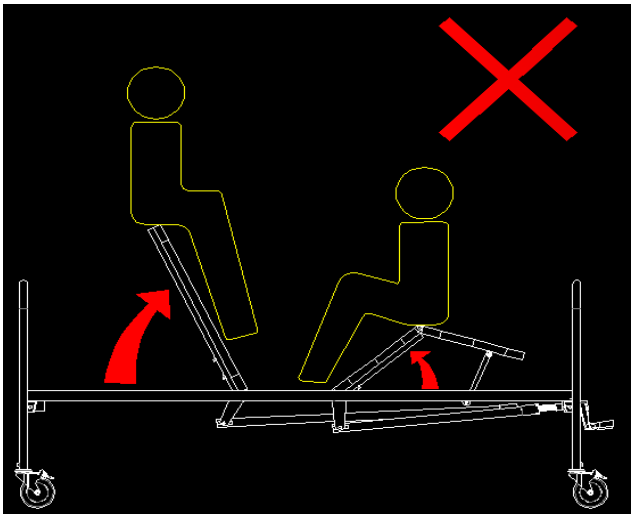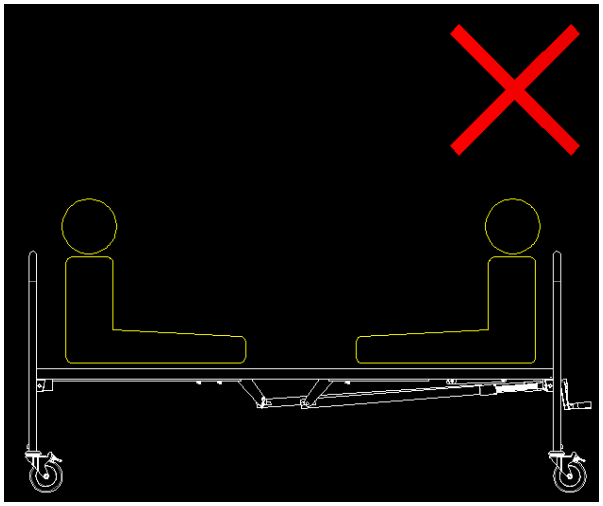10. Assembling:

- The product comes assembled ready for use.
- Remove the product from the packaging.
- Release all connections and wiring from the ties.
- Plug the power cord into a power outlet.
- Press the appropriate buttons on the control panel to adjust the backrest, leg & knee tilt as well as the desired height of the bed from the floor.
- You will need to install the rails and the head and footboard

NOTE: The assembly illustrations are indicative to make assembly easier and may differ from reality. The package includes tools for assembly.

.Head-board installation

Insert the head- board to the hole.

Turn the hook to the hole if need to lock the head-board.

Foot board installation

Insert the foot-board until heard “Click” sound.

To remove foot board, force the latch outwards and pull the foot board out.

Fifth wheel pedal

Brake pedal

Pull this lever to lay down the side rail.

Indicator Light


When enabled, it is green

When it is Green the central brake is not activated


The light turns green when the "night light" is activated.

Built-in bed controller

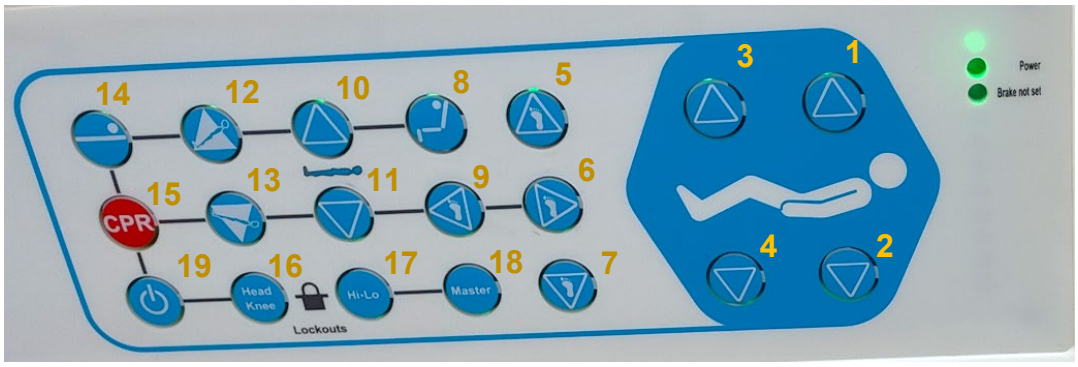1) Backrest Up
2) Backrest Down
3) Knees Up
4) Knees Down
5) Leg Lift
6) Leg Extension (Outwards)
7) Leg Fold Down
8) One-button Chair Position
9) Leg Extension Fold Down
10) Bed rising height
11) Bed lowering height
12) Reversed Trendelenburg
13) Trendelenburg
14) Electric recovery to flat position
15) Electronic CPR Button
16) Foot & Head Lock
17) Bed Height Lock
18) All Functions Lock
19) On/Off

Wired controller

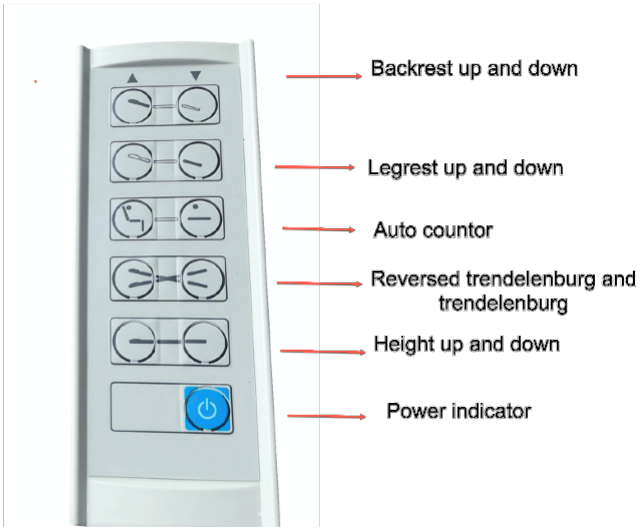Oxygen cylinder holder

Fit the oxygen cylinder holder to the headboard

It is suitable for 2,3 liters oxygen cylinders

11. Maintenance:

Check the Product after each use. The inspection is intended to investigate anything that may compromise the integrity and proper operation of the product. Check the product for any looseness, cracks, breaks, bent or damaged parts.

WARNING
In the event of a problem occurs, **DO NOT** use the product until the issue will be repaired.

12. Cleaning And Disinfection:

- You can use a wet cloth and mild soap on the surfaces that need a cleaning.
- DO NOT leave the product wet or soiled with dust.
- Do not use chemical cleaners to clean the frame. This may damage the product's coating surface.
- CAUTION**:** Do not use acidic or caustic solutions**,** Do not use solvents or stain removers.
- If you must disinfect the product, use a common, mild disinfectant, after reading the instructions provided by the disinfectant manufacturer.
- Protect the product from scratches, cuts and punctures.

13. Conditions of Storage and Transport

The product should be stored packaged under the following conditions:

- The product should be placed in a dry and ventilated place and should not be placed in an environment with high and rapid changes of temperature.
- DO NOT store the product in locations that are close to high temperature sources or direct sunlight for long periods of time.
- The product should be isolated from chemicals and corrosive substances.
- The product should be packaged to protect it from dust, and parts that may be affected or damaged by oxidation should be lightly lubricated.
- The product can be transported by any common vehicle.

14. Environmental Protection Statement:

If one day you find that your product needs to be replaced or it is no longer working for you, consider protecting the environment:
1) Do not dispose your product along with the rest of the public waste (this is also the meaning of the shown recycling sign).
2) Contact your Public Authorities and they will instruct you of the Recycling centers to which your product must be disposed.
3) Correct disposal of your product helps the protection of the environment as well as the recycling of the product’s components.
4) The electrical products due to their construction materials, if not correctly disposed, may lead to environmental and furthermore health hazards

15. Incident Reporting:

Any serious incident that occurs in relation to use of the device must be reported to the manufacturer and a competent authority of the Member State in which the user and / or patient resides.

16. Declaration Of Conformity:

We are solely responsible for declaring that the Medical Devices mentioned in this statement are of Low-Risk Class (Class I) and comply with the requirements of the European Regulation 2017/745 and where appropriate, the standards and legislation referred to.

17. Terms of Warranty:

The product is guaranteed for a period of three (3) years warranty from the purchasing date.

The warranty apply on factory defects & DOES NOT cover damages resulting from misuse, poor maintenance, modification, overuse, or non-compliance with the product's operating instructions. Also DOES NOT cover the parts that can wear out during use or over time, like the upholstery, wheels, wheel tire, safety belt, footrests or brakes.

Damages or defects caused by: natural disasters, unauthorized maintenance or repair, power supply problems (where applicable) & carriers are not covered by the warranty.

No service or spare part is covered unless it has been inspected by an authorized service department or the manufacturer of the product. Shipping costs for goods of warranty are to be paid by the user- customer or by dealer.

When the user or the provider apply on warranty request, provide a valid proof of purchase and the completed warranty form, with a photo of the product serial number and a photo or video of the problem that has occurred. In the event the user cannot provide the above documents, if we can confirm that the product or part is still within the warranty period, we may be able to honor the user's request for coverage under warranty, but we are not obligated to do so.

Warranty Form



`BUYERS DETAILS:	
FULL NAME:	
PHONE:	
PURCHASING DATE:	
INVOICE NUMBER:	
SERIAL NUMBER:	
LOT:	
DEALER DETAILS:	
FULL NAME:	
ADDRESS:	
PHONE:	
	DEALER SIGN AND SIGNATURE

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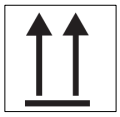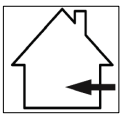Description of symbols that have been printed on the label, Master Carton and or User Manual

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