



BREEZY PREMIUM

STRONG AND DURABLE



STRONG AND DURABLE

The loyal life partner!



 B R E E Z Y
PREMIUM



Configuration with options





THE STEEL WHEELCHAIR WITH SINGULAR STRENGTH AND DURABILITY

DURABILITY

Breezy Premium, the best allied you could desire! A chair of another generation, designed in details to last overtime thanks to its incredible **rigidity, strength and ergonomic**. Its steel structure of great quality guarantees unparalleled functionality on a day to day basis, making it as a chair of total confidence.



DETAILS THAT MAKE THE DIFFERENCE

Breezy Premium offers many features to make life easier. All the elements have been taken carefully to offer a confidence chair, ultra strong and durable. To offer... a Premium product! Details as **brakes lever with ergonomic handle**, with the possibility of regulating them or the **wide unlocking lever of the footrest** make the Breezy Premium a chair which makes a strong impression.



Knee lever brake with comfortable and ergonomic handle, effective braking with minimum effort.



Detachable footrest, can swing outwards and inwards with an ergonomic, rigid and handy latch. Adjustable footplate according to wheelchair width.



THE BEST OPTION

Possibility of choosing between **3 different models** depending on the backrest you prefer (Rigid, Reclining or Half-folding). Anyway, has the option of **choose between 4 colours** plus a **wide range of accessories** to ensure the chair to meet the user's needs perfectly.



3 BACKRESTS



Rigid Backrest



Reclining Backrest



Half-folding Backrest

COLOUR RANGE

4 frame colour.



Forest Green



Ruby Red



Navy Blue



Selenium Grey

ADJUSTABLE - DISTANCE BETWEEN WHEELS

Modifying the position of the rear wheel plate we can change the distance between the front and rear wheelchair. The length of the chair can be lengthened or shortened depending on the turning radius that want to acquired. This helps to gain an more active chair, while lengthening an stable chair will be gotten.



Set forward
Std. position at the center of gravity (+25mm). More active



Set back
Center of gravity -50mm
More stable

3 SEAT HEIGHT

Breezy Premium can achieve **3 different seat height** (50,5, 48 or 45,5 cm). The height depend on the rear wheels and castors

COMFORT WITHOUT LIMITS

REDUCED TOTAL LENGTH

The unique design of the front frame and the rear wheel position give a total length of just 76.5cm without footrests. This helps the chairs to gain better access to confined spaces.



WIDE RANGE OF OPTIONS AND ACCESSORIES

Breezy Premium is adapted to the needs of different end users. It offers a wide range of options and accessories, in order to adapt the chair to any person requirements.


- Anatomical seat and backrest
- Breathable vented upholstery
- Aluminium angled adjustable footplates
- Elevating legrests
- Footboard
- Headrest
- Drum brakes
- IV pole
- Crutch holder
- Wheel lock extension
- Axle for double amputees
- Anti tips
- Transit wheels
- Oxygen tank holder
- Tray table
- Belt
- Back pack..



RIGID BACKREST

STANDARD CONFIGURATION

- Steel wheelchair with rigid backrest
- Nylon upholstery
- Folding armrest
- Standard footrest with 70° hanger, detachable and swing outwards and inwards
- Sport rear wheels
- Colours: forest green, ruby red, navy blue or anthracite grey

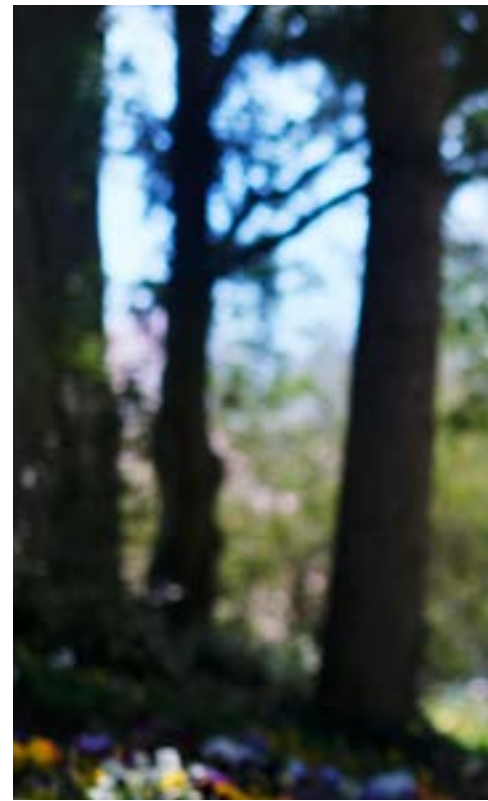

Standard configuration



**BREEZY PREMIUM WITH RIGID
BACKREST, 12" SOLID REAR WHEEL
AND BREATHABLE VENTED
UPHOLSTERY**



RECLINING BACKREST

STANDARD CONFIGURATION

- Steel wheelchair with reclining backrest till 60°
- Nylon upholstery
- T armrest, detachable in height and depth (with tools)
- Standard footrest with 70°, detachable and swing outwards and inwards.
- Sport rear wheels
- Colours: green forest, ruby red, navy blue and anthracite grey



Standard configuration with upholstery recline headrest and (optional) detachable axles



BREEZY PREMIUM OPTIONAL CONFIGURATION WITH ACCESSORIES

- 12" solid rear wheel
- T armrest, detachable in height and depth (with tools)
- Anatomical seat
- Skay headrest
- Elevating legrests



HALF-FOLDING BACKREST

STANDARD CONFIGURATION

- Steel wheelchair with half-folding backrest fo easier trasportation and storage.
- Nylon upholstery.
- Padded detachable armrest and can be folded back (with tools).
- Standard footrest with 70°, detachable and swing outwards and inwards.
- Sport rear wheels
- Colours: green forest, ruby red, navy blue and antharcite grey.



Standard Configuration



Breezy Premium
folded



**BREEZY PREMIUM WITH
HALF-FOLDING BACKREST
AND 12" SOLID REAR
WHEELS**



OPTIONS AND ACCESSORIES



UPHOLSTERY LEGREST AND FOOTPLATES



Breathable Vented upholstery



Footboard



Aluminium angel adjustable footplates



Elevating legrests

ARMRESTS



T armrests, height adjustable (with tools)



Height adjustable armrest

BACKREST



Height adjustable handles

HEADREST



Sky headrest assembled on nylon upholstery



Anatomical headrest assembled on nylon upholstery



Sky headrest assembled on anatomical headrest



Anatomical headrest assembled on anatomical backrest



ANATOMICAL OPTIONS AND BELTS



Anatomical backrest

Anatomical seat

Anatomical seat with toilet bucket

Separating wedge

Belt

Breezy Premium model with reclining backrest and anatomical seat and backrest

WHEELS



Pneumatic front castors
8" x 2"



Slid front castors
8" x 2"



Slid front castors 6"



Solid or Pneumatic rear
wheels 12"



Slid or Pneumatic spoke rear
wheels 24" or Pneumatic spoke
wheels 22"



Pneumatic mag rear wheels 24"



Quick release axle on rear wheels

BRAKES



Brake extension



Drum brake wheels available with 24" and
12" rear wheels

ACCESSORIES



Backbag



Backpack



Swing-away tray



IV Pole



Cutch holder



Oxygen bottle holder



Transit wheels



Anti-tip wheels



Double amputee kit

Other options: hemiplegia system, short & fixed footrests, antiplateled system....

TECHNICAL SPECIFICATIONS

Seat width	37 / 40 / 43 / 46 / 49 / 52 cm
Seat depth	42 cm
Seat height	50,5 cm (48 ó 45,5 cm w/22 and 6" wheels)
Backrest height	44,5 cm (82 cm reclining backrest and headrest)
Total width	Seat width + 19 cm (24" wheels)
	Seat width + 16,5 cm (12" wheels)
Total length w/24" rear wheels	103 cm (76,5 cm w/o footrest)
Total length w/12" rear wheels	95 cm (68,5 cm w/o footrest))
Total lenght reclining backrest	113 cm (87,5 cm w/o footrest)

Total lenght reclining bacrest 12" wheel	105 cm (79,5 cm w/o footrest)
Total height	94 cm (128 cm reclining model w/headrest)
Seat angle	3°
Backrest angle	8°
Colours	Green forest, ruby red, navy blue and antharcite grey
Aproved user weight	125 kg
Max. weight	From 16.5 kg (complete)*
Folded width	31 cm

* Breezy Premium Model Standard configuration, size 37.
Tolerance +/- 1 cm

For further information in the full specification, options and accessories, please refer to the order form.
All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.



Sunrise Medical S.L
Polígono Bakiola, 41
48498 Arrankudiaga
Vizcaya, España
Tel.: 902 14 24 34
Fax: 94 648 15 75
www.SunriseMedical.eu



Ver 1 - May 2018

